



SONIC FACIAL ROLLER | 4 IN 1



Skin Scan & Personal
Routine Guide



SmartSonic
Pulsation Technology



Smoothing & Tightening
Skin Massage



3D Contouring &
Firming Waves

Tighter

Skin in Minutes



**Create your personal
Skincare Routine**

and Download our GESKE App



GESKE Beauty Tech GmbH
Leipziger Platz 18
10117 Berlin Germany
info@geske.com, www.geske.com



GOOD TO KNOW

TIP 01

For best results, follow with a GESKE serum or moisturizer that is specially formulated to achieve glowing skin.

TIP 02

Always roll the device in an upward motion for best results.

TIP 03

Clean your device with a GESKE device cleanser after every use.



Smoothing & Tightening
Youthful Skin



Massage
Glowing Skin




On/Off
Press for 1-2 Seconds



Contouring & Firming
Anti-Aging




Sonic Pulsations
Soothe & Relax

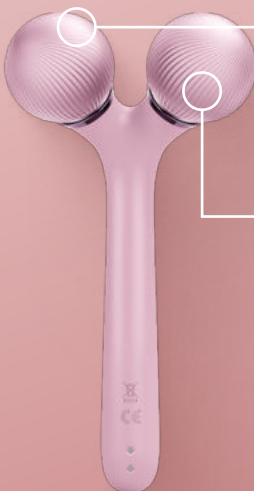

MASTER THE FIRST STEPS

Our innovative SmartAppGuided™ Sonic Facial & Body Roller | 4 in 1 is designed to visibly firm and refine the skin on your entire body, including your face.

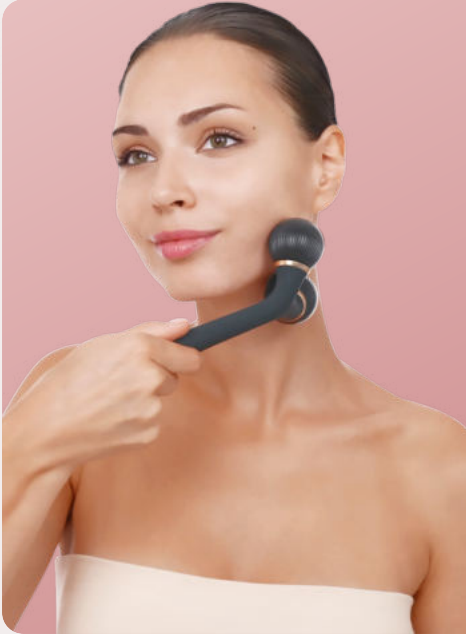

FACE & NECK MASSAGE

PREPARE

Start with a clean, make-up free face and apply your desired GESKE serum.

POWER ON


Turn on the SmartSonic Pulsation Technology by holding the power button for 2 seconds.

MESSAGE

Roll the massaging spheres in an upward motion over your décolleté, shoulders and neck. Then, move on to your face and massage your cheeks and forehead.

FINISH

Turn the device off by holding the power button for 2 seconds. Clean your device with a GESKE Device Cleanser after every use.

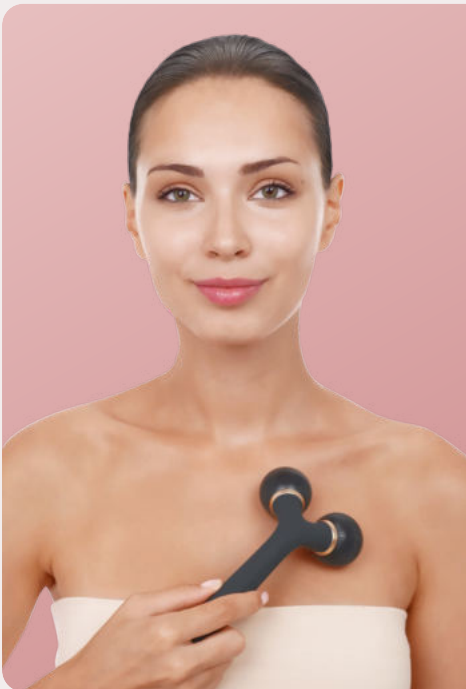

POINTS TO TARGET

NECK

Starting at the bottom of your neck, roll it in an upward motion towards the jawbone. Repeat five times on both sides of the neck. The 3D Contouring and Firming Waves enable your skin to trigger its natural renewing process.

DÉCOLLETÉ

Roll the device across the collarbone, towards the shoulder. Repeat five times on both sides of your body. The SmartSonic Pulsation Technology stimulates and lifts your skin.

CHEEKS

Roll the device from your chin to the ear, along the jawline. Repeat five times on both sides of your face. This Smoothing and Tightening Skin Massage gives you firm, pearly-smooth skin.

FOREHEAD

Starting at the eyebrows, roll the device upwards, to the hairline. Repeat five times, covering all areas of your forehead.

Create your personal Skincare Routine

and Download our GESKE App



IMPORTANT SAFETY INFORMATION. READ BEFORE USE.

Read the safety information before first use. Keep both instructions for use.

Your product at a glance

Your product is a multifunctional derma roller. It is safe to be used up to 3 times a day, please do not use more frequently. For maximal safety and ideal results, always follow your personalized Skincare Routine created by the GESKE App.

Notes on Use

Before using the device for the first time, this Safety Information and the Product Manual must be read in their entirety. Please make sure to save a copy of both as they contain important information on safe use as well as on technical and legal information requirements. The device is a skincare product for adults and not a medical device or a product with medical or other healing properties. Use of the device should always include observance of the personalized application recommendations from the GESKE Beauty App.

Specifications

Materials: ABS, silicone | **Storage temperature:** -5 - +60 °C/23 - 140 °F | **Battery:** Li-Ion 3.7 V | **Battery Capacity:** 140 mAh | **Maximum Charging Current:** 140 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 2.5 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration, 1 intensity

In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:

In case of pregnancy. The safety of use during pregnancy has not been established.

Further information on use:

Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere). This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product. Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored. If pain and/or discomfort occurs during use, discontinue use and seek medical advice. The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself. Only use the product as recommended in your personalized routine. For hygiene reasons, we recommend you do not share your product with anyone else. The product is intended for cosmetic use on the face and parts of the body, as advised in the GESKE Beauty App. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls. Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet.

Warnings

Risk of injury if used incorrectly!

Only use the product as described in the Quick Start Guide and video instructions. If you still have questions, please contact our customer support at support@geske.com.

Risk of injury due to fire/ignition!

Do not extinguish the product with water. Do not throw the product into fire.

Risk of injury due to defective batteries!

Improper handling of batteries can cause an uncontrolled exothermic reaction. In such a case, dispose of the entire product immediately and properly.

Risk of injury to children!

Keep the product out of the reach of children.

Risk of injury from electric shock!

Liquid can damage the charging circuit and endanger your health. Disconnect the product from the charging cable/circuit while cleaning.

Risk of damage from magnetic fields!

Magnetic fields can affect the mechanical and electronic components such as pacemakers, for example. Do not use the product directly over a pacemaker. Contact your treating physician for more information. Do not place cards with magnetic strips (e.g. credit cards) near the product, as the product's magnets may damage these.

Risk of damage due to changing surrounding temperatures!

Changing surrounding temperatures (e.g. after transport) may lead to the formation of condensation. In this case, do not operate the product until it has adjusted to the new surrounding temperature. Never expose your product to direct sunlight, extreme heat or boiling water.

Disconnect the charging cable from the circuit after each charging process.

Safe charging up to a surrounding temperature of 45 °C/113 °F.

Safe use up to a surrounding temperature of 50 °C/122 °F.

Risk of damage from charging too long!



Information on device disposal and WEEE. This electrical appliance was first released after August 13, 2005. In accordance with the European WEEE Directive and national regulations, the electrical appliances and electronic parts included in this delivery must not be disposed of via your household waste. At the end of the product life, they must be disposed of at a municipal collection point in order to ensure environmentally friendly disposal.

Notes on batteries / rechargeable batteries

In the context of distributing batteries and rechargeable batteries, in accordance with European regulations we are obliged to inform you, the consumer, of the following: You are legally obliged to return batteries and rechargeable batteries. After use, you can return them to us, take them to a municipal collection point, or deliver them to a local dealer free of charge. Batteries and rechargeable batteries are labelled as follows:



These symbols mean that the batteries and rechargeable batteries cannot be disposed of by you, the consumer, via household waste. Batteries and rechargeable batteries containing harmful substances are labelled with one of these symbols. The symbol consists of a crossed-out wheeled bin and the chemical symbol of the heavy metal which constitutes the main pollutant. In detail, the aforementioned symbols mean:

Pb = battery contains lead, Cd = battery contains cadmium, Hg = battery contains mercury

Notes on lithium batteries/rechargeable batteries

You can recognise lithium batteries and lithium rechargeable batteries through the chemical symbol „Li“. In addition to the information provided above, we would also like to point out the following: When disposing of lithium batteries or rechargeable batteries, either by returning them to us, taking them to a municipal collection point or handing them in to local retailers or public waste disposal companies, please ensure that they are completely discharged. We also request that you cover the contact surfaces on the battery marked „+“ and „-“ with adhesive tape to avoid the risk of fire.



GESKE Beauty Tech GmbH, Leipziger Platz 18, 10117 Berlin hereby declares that all electronic devices conform with the relevant conditions of the applicable European Directives 2014/30/EU and 2011/65/EU.

Limitation of liability

GESKE Beauty Tech GmbH and/or their retailers are liable for damages caused by intentional or grossly negligent behaviour of GESKE Beauty Tech GmbH, a representative or vicarious agent, in accordance with legal directives. With the exception of injury to life, limb or health, and breach of essential contractual obligations (cardinal obligations), GESKE Beauty Tech GmbH and/or their retailers are only liable for damages which can be attributed to intentional or grossly negligent behaviour. This also applies to indirect damages, in particular lost profits. Moreover, liability with regard to consumers due to an assumed guarantee, except in the case of intentional or grossly negligent behaviour, or where there is an injury to life, limb or health and a breach of essential contractual obligations (cardinal obligations) is limited to the typical, foreseeable damages on conclusion of the contract, and also to the level of the average damages typical for this contract. This also applies to indirect damages, in particular lost profits. With regard to companies, except in the case of injury to life, limb and health or intentional or grossly negligent behaviour by GESKE Beauty Tech GmbH, and/or its retailers liability is limited to the typical, foreseeable damages on conclusion of the contract, and also to the level of the average damages typical for this contract. This also applies to indirect damages, in particular lost profits. No liability will be accepted for damages due to: Non-compliance with this manual/improper use/Unauthorised alterations/ Technical modifications/Use of unapproved spare parts/Use of unapproved accessories Claims for liability due to the German product liability legislation remain unaffected.