



SONIC WARM & COOL MASK | 9 IN 1



Skin Scan & Personal
Routine Guide



SmartSonic
Pulsation Technology



Cryo Deep Cooling
Revitalizing Technology



Pore-Opening Deep
Warming Technology



Youthful Glow
Massage Technology



Full Facial Hyper-
Infusion Technology



Full-Spectrum
LED Light Technology



High Efficiency
Depuffing System



Eye-Nose Area
Targeted Design

2

Minute
Supercharged
Facial



**Create your personal
SkinCare Routine**

and Download our GESKE App



GESKE Beauty Tech GmbH
Leipziger Platz 18
10117 Berlin Germany
info@geske.com, www.geske.com

USE WITH GESKE DEVICES ONLY



FIRST STEPS

OPEN

Take off the magnetic outer ring and place the GESKE tissue on the device's metal surface.

ATTACH

Attach the magnetic ring to easily secure the tissue.

APPLY

Apply your favorite GESKE mask, serum or moisturizer to the tissue.



On/Off

Press for 1-2 Seconds



Massage

12 Stages



Warming

Pore-Opening



Cooling

Revitalizing




Mask Application

Super efficient



LED-Light Technology

8 Colors

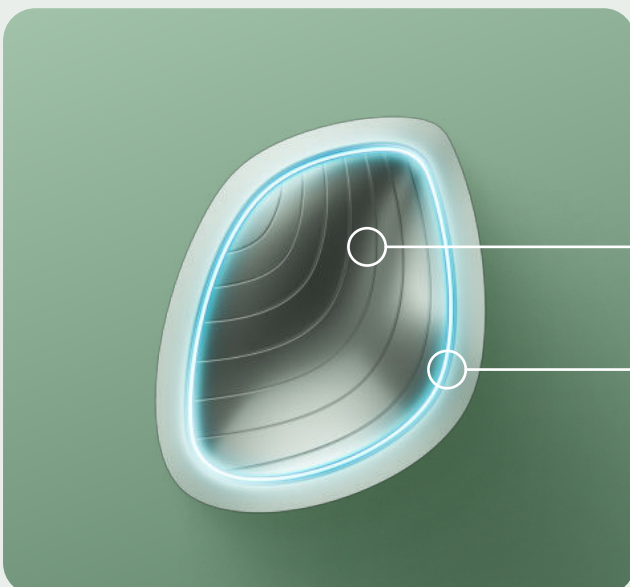

MASTER THE FIRST STEPS

Effectively apply a facial mask, fight first signs of aging or enjoy hot and cold massages – all that and much more with our SmartAppGuided™ Sonic Warm & Cool Mask | 9 in 1.



SESSION PREPARATION

PREPARE

Gently remove the ring securing the mask from the top of the device. Place the fabric sheet onto the device and secure it with the ring. Then, apply a GESKE Mask onto the fabric.

POWER ON

Turn on the SmartSonic Pulsation Technology by holding the power button for 2 seconds. Select the desired intensity by using the plus and minus buttons.

WARMING

Press the sun icon to turn on the Pore-Opening Deep Warming Technology.

LED LIGHTS

Press the power button to cycle to the desired LED color. Each color has a different nurturing effect on your skin.



2-MINUTE SUPERCHARGED FACIAL

APPLICATION

Rub the device on your face in small circles. Continue this for 90 seconds, making sure to cover all areas of the face.

EFFECT

Our proprietary Pore-Opening Deep Warming Technology helps your skin to efficiently absorb the active ingredients of your favorite GESKE mask.

COOLING

After finishing the application, press the snowflake icon to turn on the Cryo Deep Cooling Revitalizing Technology. This closes your pores, locking in moisture.

FINISH

Turn the device off by holding the power button for 2 seconds. Remove the tissue and dispose it of. Clean your device with a GESKE Device Cleanser after every usage and place the ring back on.

Create your personal SkinCare Routine

and Download our GESKE App



Advanced
App



Video
Guided



AI Powered
Skin Scan



Expert
Tips



Personal
Guide



In App
Rewards



IMPORTANT SAFETY INFORMATION. READ BEFORE USE.

Read the safety information before first use. Keep both instructions for use.

Your product at a glance

Your product is a multi-functional mask applicator. It is safe to be used up to 3 times a day, please do not use more frequently. Please do not use the Light and Pulsation mode more than once per day. For maximal safety and ideal results, always follow your personalized Skincare Routine created by the GESKE App. Download the GESKE app and get your personalized routine.

Notes on Use

Before using the device for the first time, this Safety Information and the Product Manual must be read in their entirety. Please make sure to save a copy of both as they contain important information on safe use as well as on technical and legal information requirements. The device is a skincare product for adults and not a medical device or a product with medical or other healing properties. Use of the device should always include observance of the personalized application recommendations from the GESKE Beauty App. The device is rechargeable and allows for convenient and wireless operation. It must be fully charged before the first use. All you need is the included USB charging cable. The initial charge can take up to eight hours. The cable must dock with both magnetic contacts on the device's designated metal contact points. If the charging cable does not attach by itself, connect it to the power and wait. This reactivates the magnetic effect and the contacts connect automatically. Make sure that the magnetic contacts only connect to the contact points of the device and not to other metallic surfaces.

Specifications

Materials: ABS, silicone, PC, aluminium alloy | **Storage temperature:** -5 - +60 °C/23 - 140 °F | **Battery:** Li-Ion 3.7 V | **Battery Capacity:** 800 mAh | **Maximum Charging Current:** 450 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 4 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration: 12 intensities, Heating, Cooling, LED Light Technology; 8 colors

These symbols mean that the batteries and rechargeable batteries cannot be disposed of by you, the consumer, via household waste. Batteries and rechargeable batteries containing harmful substances are labelled with one of these symbols. The symbol consists of a crossed-out wheeled bin and the chemical symbol of the heavy metal which constitutes the main pollutant. In detail, the aforementioned symbols mean:

Pb = battery contains lead, Cd = battery contains cadmium, Hg = battery contains mercury

The product is not suitable for use under the following conditions:

Do not use if you are known to be hypersensitive to heat (e.g. heat urticaria). If a feeling of (too) intense heat occurs during usage, it must be stopped immediately. The safety and effectiveness of the product has never been assessed in combination with non-GESKE skin care products and masks. We recommend using the product only with GESKE products. Use with other products at your own risk.

Further information on use:

The safety and effectiveness of the product has never been assessed in combination with non-GESKE skin care products and masks. We recommend using the product only with GESKE products. Use with other products at your own risk. The product can emit intense light, avoid applying it in the direction of or near an open eye. Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere). This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product. Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored. If pain and/or discomfort occurs during use, discontinue use and seek medical advice. The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself. Only use the product as recommended in your personalized routine. For hygiene reasons, we recommend you do not share your product with anyone else. The product is intended for cosmetic use on the face, neck and décolleté. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls. Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet.

Warnings

Risk of injury if used incorrectly!

Only use the product as described in the Quick Start Guide and video instructions. If you still have questions, please contact our customer support at support@geske.com.

Risk of injury due to fire/ignition!

Do not extinguish the product with water. Do not throw the product into fire.

Risk of injury due to defective batteries!

Improper handling of batteries can cause an uncontrolled exothermic reaction. In such a case, dispose of the entire product immediately and properly.

Risk of injury to children!

Keep the product out of the reach of children.

Risk of injury from electric shock!

Liquid can damage the charging circuit and endanger your health. Disconnect the product from the charging cable/circuit while cleaning.

Risk of damage from magnetic fields!

Magnetic fields can affect the mechanical and electronic components such as pacemakers, for example. Do not use the product directly over a pacemaker. Contact your treating physician for more information. Do not place cards with magnetic strips (e.g. credit cards) near the product, as the product's magnets may damage these.

Risk of damage due to changing surrounding temperatures!

Changing surrounding temperatures (e.g. after transport) may lead to the formation of condensation. In this case, do not operate the product until it has adjusted to the new surrounding temperature. Never expose your product to direct sunlight, extreme heat or boiling water.

Safe use up to a surrounding temperature of 50 °C/122 °F.

Risk of damage from charging too long!

Disconnect the charging cable from the circuit after each charging process.

Safe charging up to a surrounding temperature of 45 °C/113 °F.



Information on device disposal and WEEE. This electrical appliance was first released after August 13, 2005. In accordance with the European WEEE Directive and national regulations, the electrical appliances and electronic parts included in this delivery must not be disposed of via your household waste. At the end of the product life, they must be disposed of at a municipal collection point in order to ensure environmentally friendly disposal.

Notes on batteries / rechargeable batteries

In the context of distributing batteries and rechargeable batteries, in accordance with European regulations we are obliged to inform you, the consumer, of the following: You are legally obliged to return batteries and rechargeable batteries. After use, you can return them to us, take them to a municipal collection point, or deliver them to a local dealer free of charge. Batteries and rechargeable batteries are labelled as follows:



These symbols mean that the batteries and rechargeable batteries cannot be disposed of by you, the consumer, via household waste. Batteries and rechargeable batteries containing harmful substances are labelled with one of these symbols. The symbol consists of a crossed-out wheeled bin and the chemical symbol of the heavy metal which constitutes the main pollutant. In detail, the aforementioned symbols mean:

Pb = battery contains lead, Cd = battery contains cadmium, Hg = battery contains mercury

Notes on lithium batteries/rechargeable batteries

You can recognise lithium batteries and lithium rechargeable batteries through the chemical symbol „Li“. In addition to the information provided above, we would also like to point out the following: When disposing of lithium batteries or rechargeable batteries, either by returning them to us, taking them to a municipal collection point or handing them in to local retailers or public waste disposal companies, please ensure that they are completely discharged. We also request that you cover the contact surfaces on the battery marked „+“ and „-“ with adhesive tape to avoid the risk of fire.



GESKE Beauty Tech GmbH, Leipziger Platz 18, 10117 Berlin hereby declares that all electronic devices conform with the relevant conditions of the applicable European Directives 2014/30/EU and 2011/65/EU.

Limitation of liability

GESKE Beauty Tech GmbH and/or their retailers are liable for damages caused by intentional or grossly negligent behaviour of GESKE Beauty Tech GmbH, a representative or vicarious agent, in accordance with legal directives. With the exception of injury to life, limb or health, and breach of essential contractual obligations (cardinal obligations), GESKE Beauty Tech GmbH and/or their retailers are only liable for damages which can be attributed to intentional or grossly negligent behaviour. This also applies to indirect damages, in particular lost profits. Moreover, liability with regard to consumers due to an assumed guarantee, except in the case of intentional or grossly negligent behaviour, or where there is an injury to life, limb or health and a breach of essential contractual obligations (cardinal obligations) is limited to the typical, foreseeable damages on conclusion of the contract, and also to the level of the average damages typical for this contract. This also applies to indirect damages, in particular lost profits. With regard to companies, except in the case of injury to life, limb and health or intentional or grossly negligent behaviour by GESKE Beauty Tech GmbH, and/or its retailers liability is limited to the typical, foreseeable damages on conclusion of the contract, and also to the level of the average damages typical for this contract. This also applies to indirect damages, in particular lost profits. No liability will be accepted for damages due to: Non-compliance with this manual/Improper use/Unauthorised alterations/ Technical modifications/Use of unapproved spare parts/Use of unapproved accessories Claims for liability due to the German product liability legislation remain unaffected.